



Quick Start Guide

USA

neuronode@controlbionics.com
+1 (855) 831-7521
745 Center Street, Suite 303
Milford OH, 45150

AUSTRALIA

neuronode@controlbionics.com
+61 3 9897 3576
171 Union Road, Suite G.03
Surrey Hills VIC 3127


WHATS IN THE BOX



Cosmos Connect x1



Micro USB Cable x1

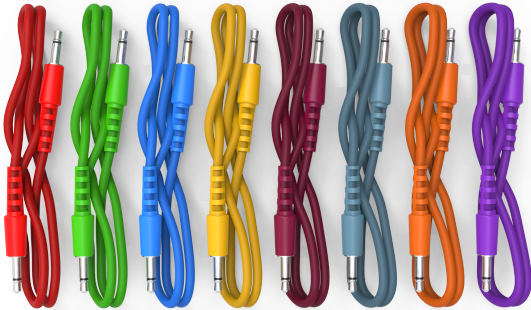

3.5mm Colour Coded Cables x 8

COSMOS CONNECT SETUP




Connect the micro USB cable into the Cosmos Connect unit and then plug the other end into any available USB slot in your Windows device.


DOWNLOAD

On your Windows device open your preferred browser, type in the link below and download the latest Cosmos Connect software. The software will be downloaded in a ZIP file to your general downloads folder.

www.controlbionics.com/cosmos-connect-download/

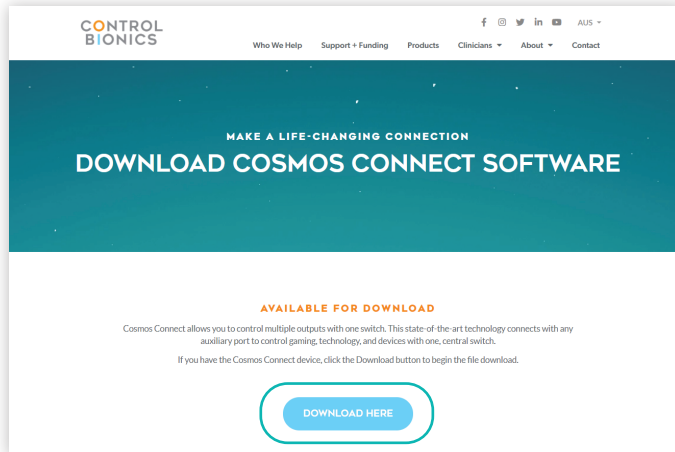


INSTALL

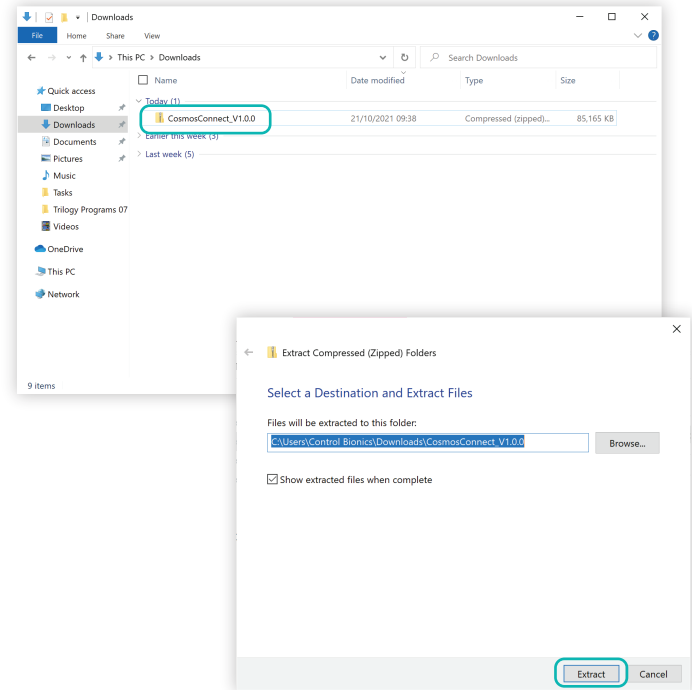


Step 1

To install the Cosmos Connect software, open the downloads folder and right click or hold the CosmosConnect ZIP file and select "Extract All".

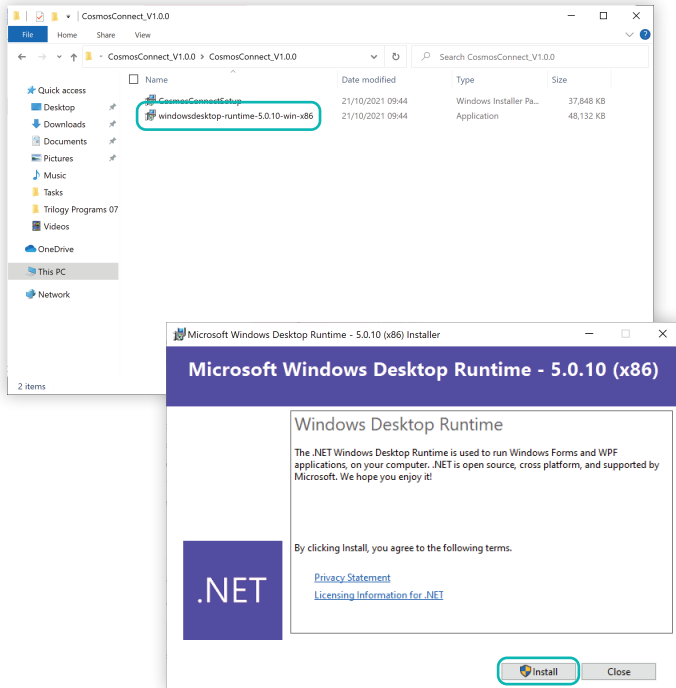


Step 2

Open the CosmosConnect folder. Double click "Windowsdesktop-runtime" and then click Install.

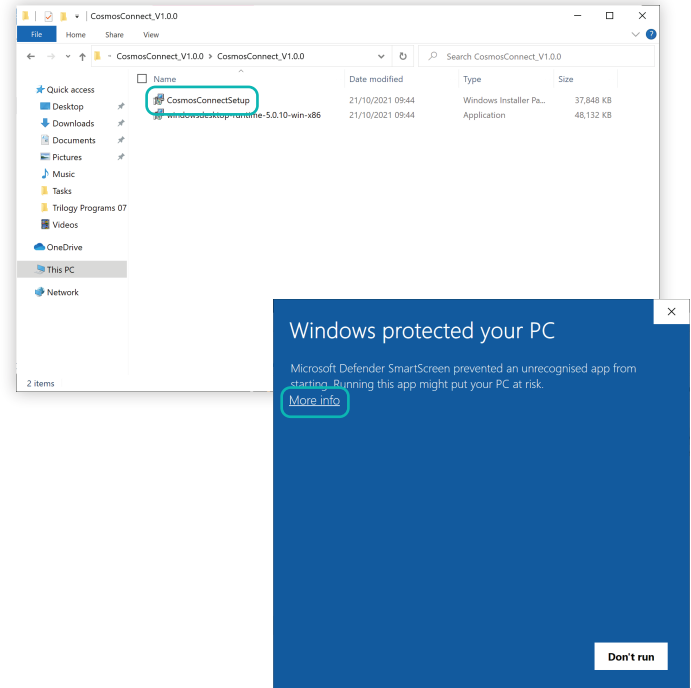

Step 3

Next Double click the Cosmos Connect Setup application. On the popup click "More Info".

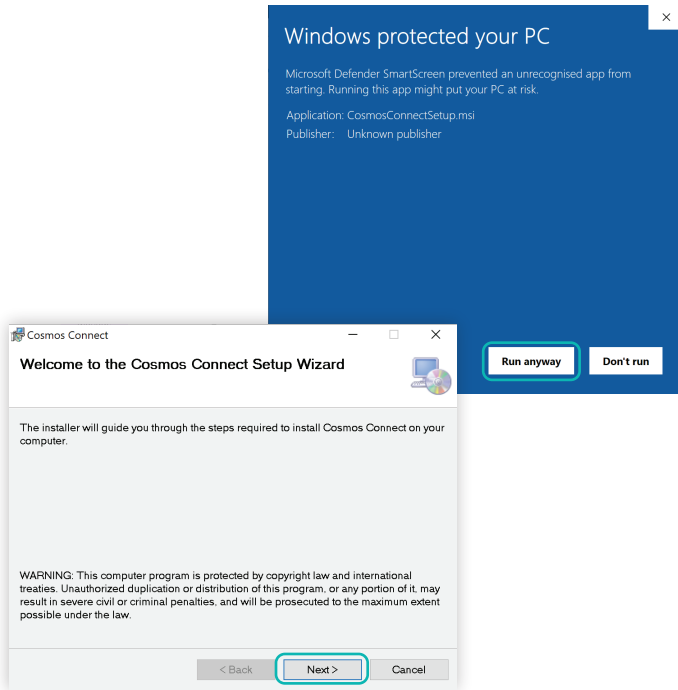


Step 4

Click "Run Anyway", Then follow the Cosmos Connect Setup Prompts. Once complete, close the setup wizard and a Cosmos Connect desktop icon should be installed.

LAUNCH SOFTWARE



To launch Cosmos Connect double click this icon located on your Trilogy's desktop.

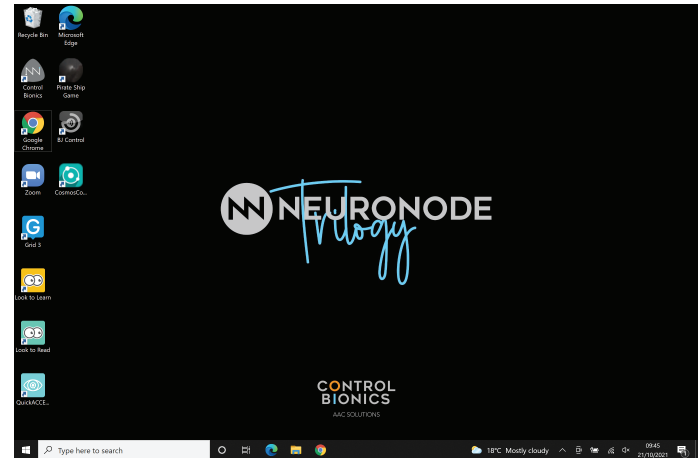

SELECTING YOUR APPS

Now that your Cosmos Connect is all setup, begin by pressing "Start" and then follow the prompts until you arrive at the "Choose Presets" page.

Here you can choose what apps you intend to use.



Apple iOS

Control and interact with your Apple iOS device, such as a iPhone or iPad.

Apple name and logo are property of Apple Inc.



Android OS

Control and interact with your Android OS device, such as a Galaxy Phone or Galaxy Tab.

Android name and logo are property of Google LLC.

Use of these names and logos does not imply affiliation with or endorsement by their respective owners.



XBOX Adaptive Controller

Control and play a variety of Xbox games. You will need to choose the specific Xbox games available to you from the pre-set menu.

Xbox name and gaming graphics are property of Microsoft Inc.



Obi Robot

Control an Obi Feeding Robot to independently select your food and feed yourself.



Switch Toys

Control and play with your 1, 2 or 3 switch AAC adapted toys.

COSMOS CONNECT DESKTOP

Menu

Clicking this icon will open the settings menu.

Access method indicator

This displays the current access method that is being used on this specific page.

Close

Clicking this X will exit the Cosmos software.



Help

Opens FAQ and support information.

Switch output indicator

Flashes green to indicate that a action have been triggered on the corresponding port.

Page timeout

A countdown timer that automatically sends you back a page once the timer has run out.

Play / pause

Stops the scanning of cells and the page timeout.

Wiring Diagram

Indicates how the the colour coded cables need to be connected into the Cosmos and output device.

Home

Sends you back to your home page.

Notes

Opens a text panel where you can leave setup notes for individual pages.

MENU

Access Method

Here you can change and adjust your method of interacting with buttons for a particular page.

Page settings

Here you can adjust the layout and features of a particular page such as changing the number of buttons or increasing the time before page timeout.

Button Design

Here you can directly change how an individual button functions such as changing the button image.

(This features is intended for advanced users only)

Change Preset

This will send you back to the change presets page where you can add or remove apps from your Cosmos Connect home page.

Advanced

Here you can change advanced features.

About

Additional information about the Cosmos Connect and Control Bionics.

